


Sleep and positioning

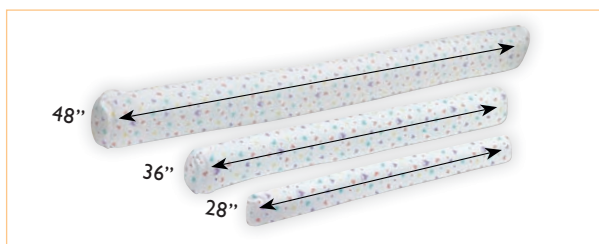

Bendy Bumper Family

The bendable positioning aid

The Bendy Bumper and the B³ Better Bendy Bumper are bendable positioning aids that stay in place until reshaped by the caregiver. They can be molded into the desired shape by the caregiver to promote containment and to provide a reflex stimulus for extremity extension and subsequent flexion recoil. This unique molding capability is due to a safe, hermetically-sealed lead rod in the center of a thick, spongy foam body.

For the Baby Bendy, the Bendy and the Super Bendy, the hermetically-sealed rod can be removed from the foam body enabling the rod to be cleaned and reused for multiple patients. The foam body is single patient use only. Two types of covers are available: washable cloth covers for multiple patient use and disposable covers for single patient use.

The Bendy Bumper is available in three sizes:

- The Baby Bendy is specially designed for infants < 1200 grams.
- The Bendy Bumper is designed for infants > 1200 grams and is appropriate for all preterm and term infants.
- The Super Bendy is designed for infants > 1800 grams. The sturdier rod is perfect for positioning multiples, larger infants and pediatric patients.

The B³ Better Bendy Bumper is a completely reusable version of the Bendy Bumpers described above.

Designed with the same patented lead rod structure, the B³ is enclosed with a permanent, wipeable urethane (Dermaplush) cover. The washable urethane cover makes the entire B³ reusable for multiple patients.

The B³ Better Bendy Bumper can be wiped down in accordance to each individual hospital unit's cleaning protocol and can be reused for multiple patients. The cloth covers can be laundered and reused; the disposable covers are designed for single patient use.

The B³ Better Bendy Bumper is available in two sizes – Better Bendy Bumper 28 inch and Better Bendy Bumper 36 inch.

Baby Bendy

Item #	Description	Qty.
92065-B	Complete unit with rod	6/case
92065-B1	Replacement foam body	6/case
92065-B2	Replacement cloth covers	12/case
92065-B12	Foam body with cloth cover	6/case
92065-BR	Baby bendy rod	6/case

Bendy Bumper

Item #	Description	Qty.
92065-A	Complete unit with rod	6/case
92065-A1	Replacement foam body	6/case
92065-A2	Replacement cloth covers	12/case
92065-A3	Disposable covers	50/case
92065-A12	Foam body with cloth cover	6/case
92065-AR	Bendy bumper rod	6/case

Super Bendy

Item #	Description	Qty.
92065-C	Complete unit with rod	6/case
92065-C1	Replacement foam body	6/case
92065-C2	Replacement cloth covers	12/case
92065-C3	Disposable covers	50/case
92065-C12	Foam body with cloth cover	6/case
92065-CR	Bendy bumper rod	6/case

B³ Better Bendy Bumper

Item #	Description	Qty.
B³ extra small – 28" long		
1029244	Complete unit with washable cover	6/case
92065-B2	Replacement cloth covers	12/case
92065-B3	Disposable covers	50/case
B³ small – 36" long		
1029243	Complete unit with washable cover	6/case
92065-A2	Replacement cloth covers	12/case
92065-A3	Disposable covers	50/case